

EXHIBIT A

酵 素

Enzyme for Biochemistry

Reagents


MERCK

**ENZYMES
for**

Biochemical research



Gene research

Food analysis

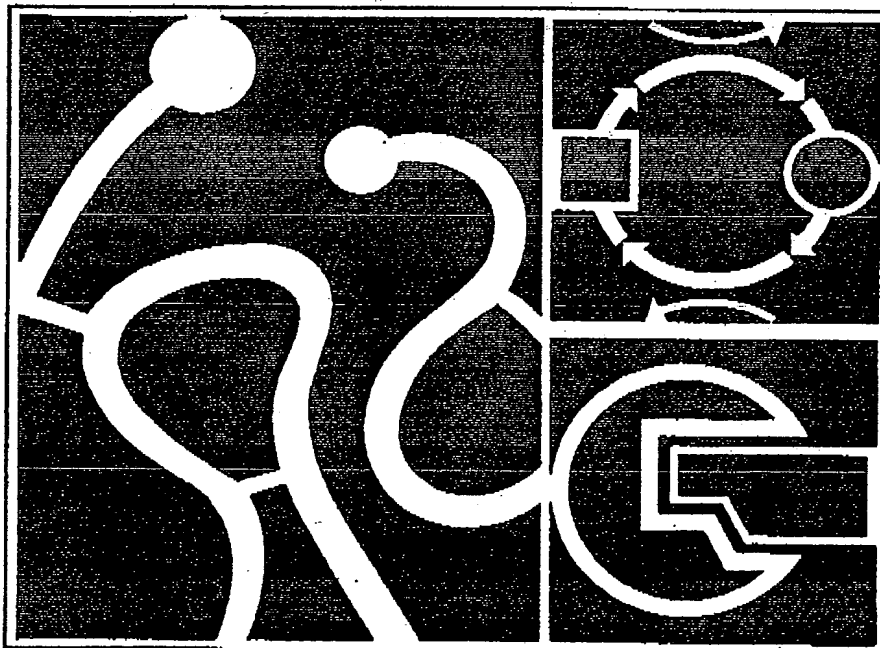

1993年 9 月

MERCK

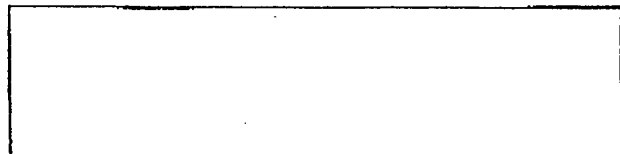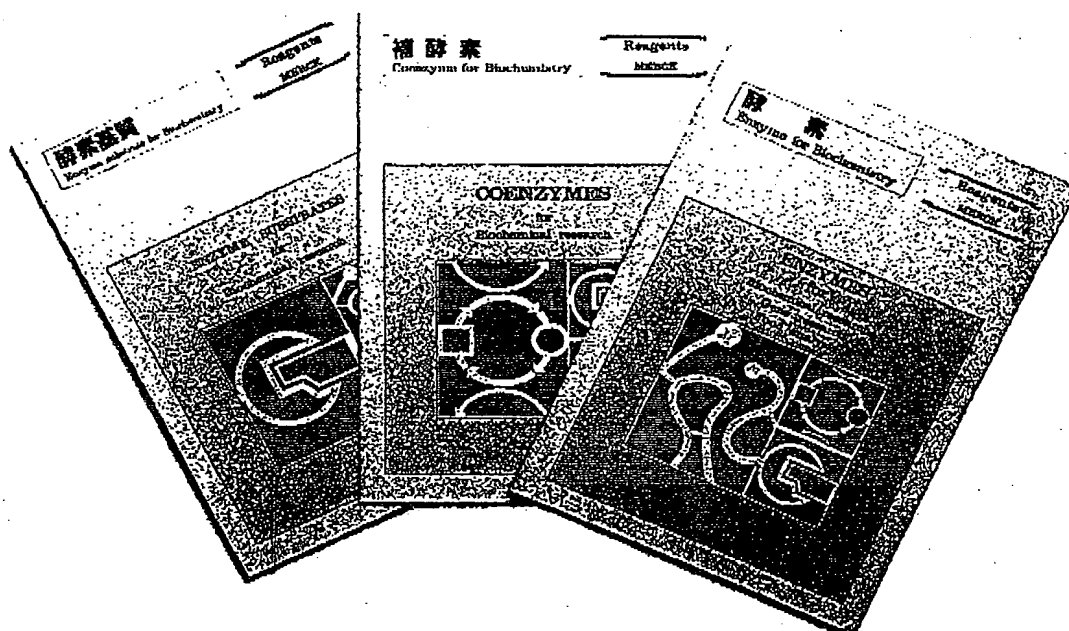

関東化学株式会社
試薬事業本部

〒103 東京都中央区日本橋本町3-2-8 03(3663)7631
〒541 大阪市中央区瓦町2-5-1 06(222)2796
〒812 福岡市博多区山王1-1-32 博多塔池ビル 092(416)9361

24686 Ribonuclease A (from bovine pancreas)

lyophilised

25 U/mg for biochemistry

EC 3.1.27.5

リボヌクレアーゼA (別名: RNase A)

反応 リボ核酸を加水分解し、2', 3'-環状ピリミジンヌクレオチドを経由し、3'-ピリミジンヌクレオチドを生成

起源 ウシ膵臓

形状 凍結乾燥品

活性 25 U/mg 以上 *no less than 25 U/mg*

Test conditions

1.4 ml Cytidine-2', 3'-cyclophosphate 22.2 mmol/l (resp. cytidine-2', 3'-cyclophosphate barium salt 11.1 mmol/l), dissolved in NaCl 0.1 mol/l (EDTA 0.1 mmol/l)

0.1 ml Ribonuclease (40 µg/ml), dissolved in NaCl 0.1 mol/l (EDTA 0.1 mmol/l)

Temperature: 25°C. Adjust to pH 7.1 and time 2 minutes. Titration of the liberated phosphoric acid groups with NaOH 0.01 mol/l by means of a pH-stat at pH 7.1.

1 U catalyses the formation of 1 µmole of phosphoric acid groups per minute under test conditions.

共存酵素活性

Deoxyribonuclease not detectable

保存 0°C - +6°C

安定性 ディープフリーザー中で、-20°Cで乾燥保存した場合、12ヶ月以内に著しい活性の低下は認められない

包装・価格 100 mg 13,600

7686 Saccharase (from yeast)

lyophilised 300000 U/vial for biochemistry

EC 3.2.1.26

サッカラーゼ

(インベルターゼ, β-フルクトシダーゼ)

系統名 β-D-Fructofuranoside fructohydrolase

反応 スクロースなどのβ-D-フルクトフラノシドの非還元性のβ-D-フルクトフラノシド残基末端を加水分解

起源 酵母

形状 凍結乾燥品

活性 約 200 U/mg

Test conditions

8.9 ml Acetate buffer 0.1 mol/l, pH 4.5

1.0 ml Sucrose 1 mol/l, dissolved in redistilled water

0.1 ml Saccharase (1 mg/ml), dissolved in redistilled water

Incubate exactly 3 resp. 6 minutes at 25°C. To stop the reaction, add 0.2 ml of the reaction mixture to 2.0 ml Tris buffer 0.1 mol/l, dissolved in redistilled water.

For the determination of the liberated glucose, mix

2.0 ml Determination reagent (phosphate buffer 0.12 mol/l, pH 7.0; NaCl 0.15 mol/l; NAD 1.1 mmol/l; 5 U/ml glucose dehydrogenase; 0.1 U/ml mutarotase)

0.2 ml Reaction mixture (stopped with Tris buffer)

Temperature: 25°C. Measurement of ΔA at 366 nm at the end of the determination reaction (ca. 10-15 minutes) against blank.

Extinction coefficient (NADH): $\epsilon_{366} = 3.4 \text{ cm}^2/\mu\text{mole}$.

1 U catalyses the formation of 1 µmole of glucose per minute under test conditions.

保存 0°C - +6°C

安定性 4°Cで乾燥保存した場合、12ヶ月以内に著しい活性の低下は認められない

包装・価格 1 pack 16,400

EXHIBIT B


BOEHRINGER MANNHEIM

Biochemica

Biochemicals

Catalog 1996 '97

研究用試薬カタログ

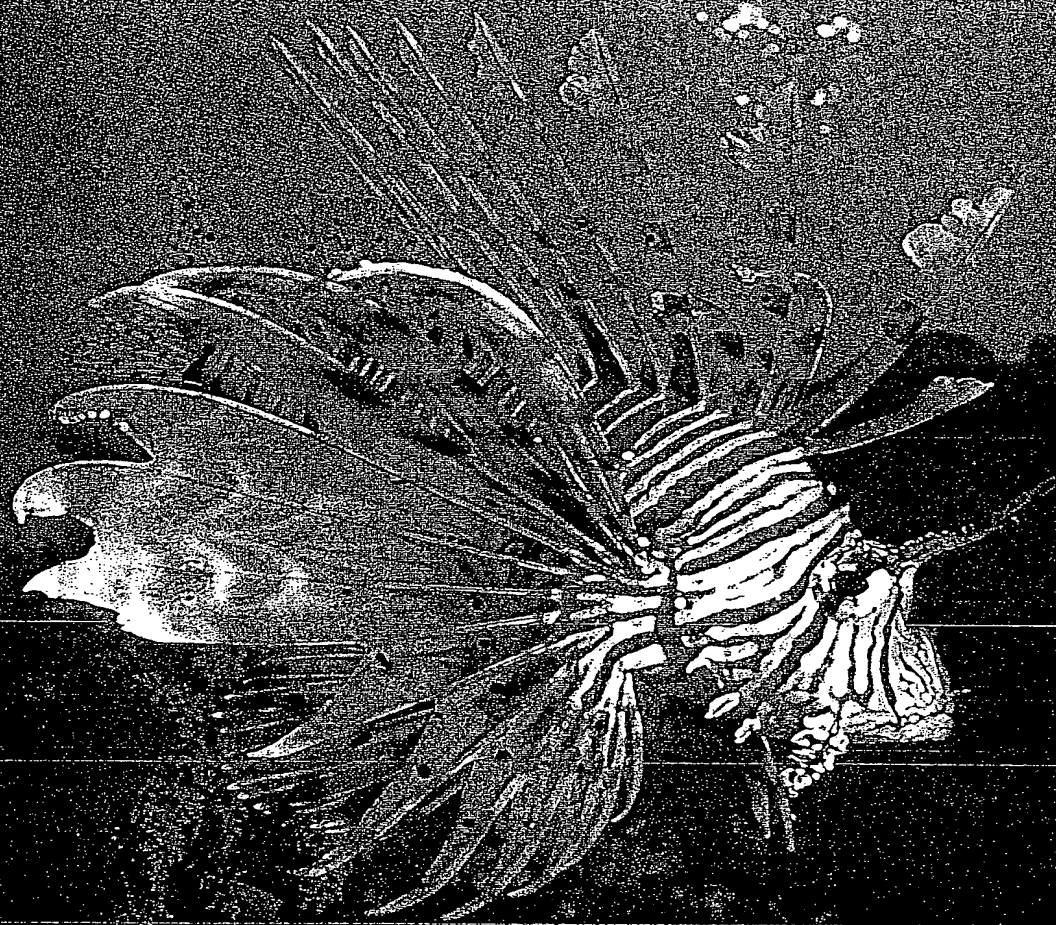

切断とシーケンス用プロテアーゼ

Proteases for Cleavage and Sequencing

プロテアーゼ

製 品 名	アプリケーション	特 異 性
トロンビン Thrombin 抽出源：ヒト血漿 凍結乾燥 緩衝剤で安定化、pH6.9 EC 3.4.21.5	血液凝固、医学研究。蛋白質構造研究及び生化学研究。	Argのカルボキシル基側でペプチド、エステル結合を特異的に切断するセリンエンドペプチダーゼ。
トリプシン Trypsin 抽出源：パン臓器 EC 3.4.21.4	蛋白質の分解と組織の分散。	塩基性アミノ酸ArgとLysのカルボキシル基側で蛋白質とペプチドを特異的に加水分解するセリンエンドペプチダーゼ。

切断とクェンス用プロテアーゼ
Proteases for Cleavage and Sequencing
プロテアーゼ

特 徴	阻 害 剤	製品番号	包装単位	希望価格
比活性：約120U/mg 酵素タンパク(Chromozym® THを基質として、25℃で測定)IU≒6.3NIH-unitsに相当 共存酵素活性：Factor Xa<3% 分子量：約 33.6kD 至適pH：8.2-9.0	DFP, TLCK, PMSF ベンザミジン, α ₁ - アンチトリプシン, α ₂ -マクログロブリン, アンチトロンピン, III-ヘパリン, ヒ ルジン, APMSE	602 400	20 U	¥10,000
形状：凍結乾燥、結晶化トリプシンより調製、塩類は含まない 比活性：約110U/mg凍結乾燥品(Chromozym® TRY ⁺ を基質として25℃で測定)≒約40U/mg凍結乾燥品(ベンゾイル-L-アルギニンエチルエステルを基質として25℃で測定) 分子量：23.5kD 至適pH：8.0	DFP, TLCK, PMSF ロイペプチン、大豆 トリプシン、インヒビ ター、卵白トリプシ ン、インヒビター、 α ₂ -マクログロブリン α ₁ -アンチトリプシ ン, APMSE、アンチ パイン	109-819 109-827	500 mg 2 g	¥ 5,000 ¥17,400

about 110 U/mg freeze-dry product (measured at 25°C using
Chromozym® TRY⁺ as a substrate) ÷ about 40 U/mg
freeze-dry product (measured at 25°C using
benzoyl-L-arginine ethyl ester as a substrate)

EXHIBIT C



MP Biomedicals, LLC.
Formerly ICN Biomedicals, Inc.

1263 S. Chillicothe Rd.
Aurora, Ohio 44202

Telephone: 330/562-1500
Toll Free: 800/854-0530
Fax: 330/562-1987
mailto:blotech@mpbio.com
web: www.mpbio.com


A Member of the American Society For Quality

ImmunO™

Nuclease(r)
Catalog #: 32035
Lot #: Typical

The *Staphylococcal* derived nuclease cleaves the 5'-phosphoryl ester bond of nucleic acids.

- Source:** *Staphylococcus aureus* gene expressed in recombinant *E. coli*
- Form:** Lyophilized in vials of 5 mg or 50 mg. The enzyme is also available in mega-unit quantities as a frozen liquid.
- Purity:** >99%, protein basis. Single band on SDS gradient gel electrophoresis on staining with Coomassie Blue R-250.
- Specific Activity:** 7,000 - 10,000 Units/mg (generally 8,000 Units/mg). One unit is equivalent to a change in A_{260} of 1.0 after 30 minutes at pH 8.8 and 37°C of a reaction mixture of acid soluble polynucleotides from native NDA. One μ molar unit = 85 A_{260} units.
- Properties:** The optimum pH for both RNase and DNase activity is between 9.0 and 10.0, and is dependent on the concentration of calcium ions. At high pH, less Ca^{2+} is required. The presence of Sr^{2+} results in high DNase activity and loss of RNase activity(1)
- Assay:** The assay of *Staphylococcus* nuclease is based on an increase in absorbance at 260 nm which accompanies the hydrolysis of nucleic acids⁽²⁾. This change has been found to correlate well with other changes reflecting DNA hydrolysis, such as an increase in viscosity, number of secondary phosphate groups which are liberated, and the production of acid soluble nucleotides⁽²⁾.
- Stability:** The enzyme is particularly stable; Stability extends to pH values as low as 0.1. At a concentration of 15 mg/ml there is no significant loss of activity after 20 minutes at 100°C.


Approved by: Joseph Dietz, Ph.D.
Quality Control Director

Control #

EXHIBIT D



Enzyme string search

Search results

Search-string [ribonuclease A]

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.4.- Phosphoric diester hydrolases.

3.1.4.22 - Transferred entry: 3.1.27.5.

PROTEIN NAME (2rsa): Ribonuclease a

PROTEIN NAME (3rsa): Ribonuclease a

PROTEIN NAME (4rsa): Ribonuclease a (joint neutron and x-ray)

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.26.- Endoribonucleases producing 5'-phosphomonoesters.

3.1.26.2 - Ribonuclease alpha.

DESC: Ribonuclease alpha.

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.27.- Endoribonucleases producing other than 5'-phosphomonoesters.

3.1.27.5 - Pancreatic ribonuclease.

PROTEIN NAME (1a2w): Crystal structure of a 3d domain-swapped dimer of bovine pancreatic ribonuclease a

PROTEIN NAME (1a5p): C[40,95]a variant of bovine pancreatic ribonuclease a

PROTEIN NAME (1a5q): P93a variant of bovine pancreatic ribonuclease a

PROTEIN NAME (1afk): Crystal structure of ribonuclease a in complex with 5'-diphosphoadenosine-3'-phosphate

PROTEIN NAME (1afi): Ribonuclease a in complex with 5'-diphosphoadenosine 2'-phosphate at 1.7 angstrom resolution

... etc

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.27.- Endoribonucleases producing other than 5'-phosphomonoesters.

3.1.27.10 - rRNA endonuclease.

PROTEIN NAME (1de3): Solution structure of the cytotoxic ribonuclease alpha-sarcin

PROTEIN NAME (1r4y): Solution structure of the deletion mutant delta(7-22) of the cytotoxic ribonuclease alpha-sarcin

Number of entries matching your search string: 4

Total number of EC entries scanned: 4327



Enzyme string search

Search results

Search-string [ribonuclease T1]

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.27.- Endoribonucleases producing other than 5'-phosphomonoesters.

3.1.27.3 - Ribonuclease T(1).

OTHER NAME(S): *Ribonuclease T1.*

PROTEIN NAME (1bir): *Ribonuclease t1, phe 100 to ala mutant complexed with 2' gmp*

PROTEIN NAME (1bvi): *Ribonuclease t1 (wildtype) complexed with 2'gmp*

PROTEIN NAME (1det): *Ribonuclease t1 carboxymethylated at glu 58 in complex with 2'gmp*

PROTEIN NAME (1fys): *Ribonuclease t1 v16c mutant*

... etc

Number of entries matching your search string: 1

Total number of EC entries scanned: 4327

Enzyme string search





Enzyme string search

Search results

Search-string [ribonuclease T2]

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.27.- Endoribonucleases producing other than 5'-phosphomonoesters.

3.1.27.1 - Ribonuclease T(2).

OTHER NAME(S): Ribonuclease T2.

Number of entries matching your search string: 1

Total number of EC entries scanned: 4327

Enzyme string search





Enzyme string search

Search results

Search-string [ribonuclease U2]

3.-.-.- Hydrolases.

3.1.-.- Acting on ester bonds.

3.1.27.- Endoribonucleases producing other than 5'-phosphomonoesters.

3.1.27.4 - Ribonuclease U(2).

OTHER NAME(S): *Ribonuclease U2.*

PROTEIN NAME (1rtu): *Ustilago sphaerogena ribonuclease u2*

Number of entries matching your search string: 1

Total number of EC entries scanned: 4327

Enzyme string search





Enzyme string search

Search results

Search-string [phosphodiesterase I]

3.-.-.- Hydrolases.
3.1.-.- Acting on ester bonds.
3.1.4.- Phosphoric diester hydrolases.
3.1.4.1 - Phosphodiesterase I.

DESC: Phosphodiesterase I

3.-.-.- Hydrolases.
3.1.-.- Acting on ester bonds.
3.1.4.- Phosphoric diester hydrolases.
3.1.4.3 - Phospholipase C.

OTHER NAME(S): Lipophosphodiesterase I.

3.-.-.- Hydrolases.
3.1.-.- Acting on ester bonds.
3.1.4.- Phosphoric diester hydrolases.
3.1.4.4 - Phospholipase D.

OTHER NAME(S): Lipophosphodiesterase II.

3.-.-.- Hydrolases.
3.1.-.- Acting on ester bonds.
3.1.4.- Phosphoric diester hydrolases.
3.1.4.17 - 3',5'-cyclic-nucleotide phosphodiesterase.

PROTEIN NAME (1ip2): Molecular docking of competitive **phosphodiesterase** inhibitor, 4-[3-(cyclopentyloxy)-4-methoxyphenyl]-2-pyrrolidinone, rolipram

Number of entries matching your search string: **4**

Total number of EC entries scanned: **4327**

Enzyme string search





Enzyme string search

Search results

Search-string [nuclease P1]

3.-.-.- *Hydrolases.*

3.1.-.- *Acting on ester bonds.*

3.1.30.- *Endoribonucleases active with either ribo- or deoxyribonucleic*

3.1.30.1 *Aspergillus nuclease S(1).*

COMMENTS: *Penicillium citrinum nuclease P1.*

Number of entries matching your search string: 1

Total number of EC entries scanned: 4327

Enzyme string search





Enzyme string search

Search results

Search-string [nuclease S1]

3.-.-.- *Hydrolases.*

3.1.-.- *Acting on ester bonds.*

3.1.30.- *Endoribonucleases active with either ribo- or deoxyribonucleic*

3.1.30.1 *Aspergillus nuclease S(1).*

OTHER NAME(S): *Aspergillus nuclease S1.*

OTHER NAME(S): *Endonuclease S1.*

OTHER NAME(S): *Deoxyribonuclease S1.*

Number of entries matching your search string: 1

Total number of EC entries scanned: 4327

Enzyme string search



**This Page is Inserted by IFW Indexing and Scanning
Operations and is not part of the Official Record**

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

- ☐ **BLACK BORDERS**
- ☐ **IMAGE CUT OFF AT TOP, BOTTOM OR SIDES**
- ☐ **FADED TEXT OR DRAWING**
- ☐ **BLURRED OR ILLEGIBLE TEXT OR DRAWING**
- ☐ **SKEWED/SLANTED IMAGES**
- ☐ **COLOR OR BLACK AND WHITE PHOTOGRAPHS**
- ☐ **GRAY SCALE DOCUMENTS**
- ☐ **LINES OR MARKS ON ORIGINAL DOCUMENT**
- ☒ **REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY**
- ☐ **OTHER:** _____

IMAGES ARE BEST AVAILABLE COPY.

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.